

Blue Tip™ (Cushion Soft Tip)

Endotracheal Tube with cushion soft tip

(Nasal/Oral)


Cuffed | SMD 772

- Endotracheal Tube - Tip Biru
- Tubo endotraqueal - Punta azul
- Tubo Endotraqueal - Ponta Azul
- Tube endotrachéal - Embout bleu

Conforming to ISO 5361



PHthalate FREE

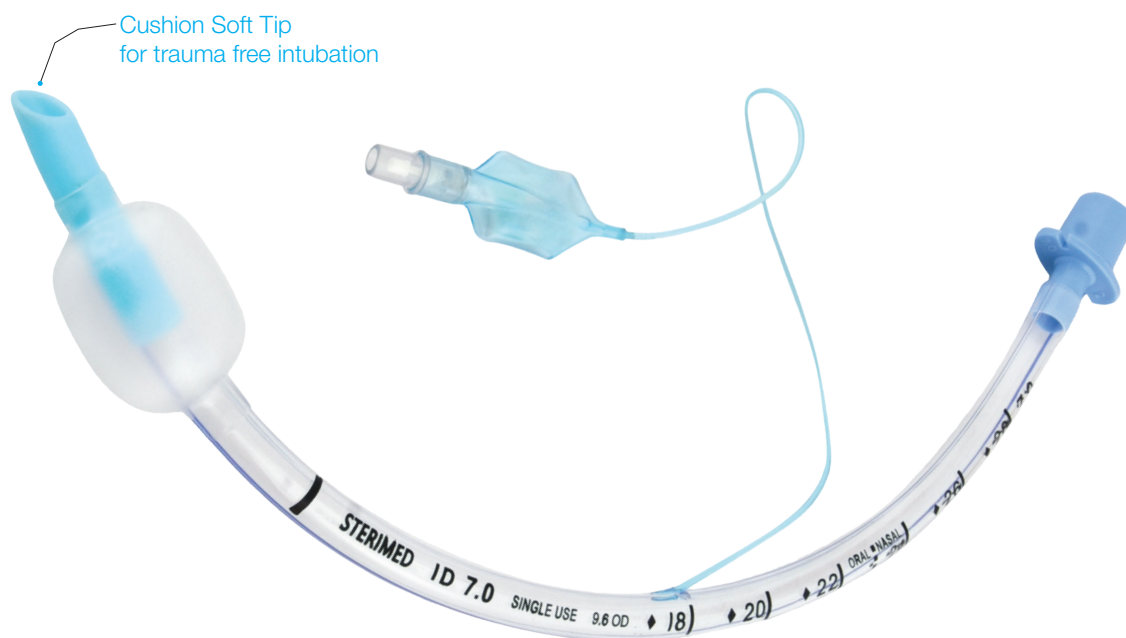

- The Tip made from a special thermoplastic elastomer, the tube offers flexibility and durability, ensuring safe and effective intubation procedures.
- Specifically designed for nasal intubation, the tube's cushion soft tip minimizes the risk of injury, prioritizing patient safety and comfort.
- Used to maintain a patent airway and mechanically ventilate patients
- Tube made from non-toxic, non-irritant medical grade clear PVC
- Kink resistant, thin-walled, thermo-sensitive tube, softens at body temperature
- Smooth inside/outside surface for ease in intubation, extubation and suctioning
- Smooth polished murphy eye, distal end and tapered balloon end for atraumatic intubation and extubation
- Kink resistant inflation tube to ensure patient safety during cuff inflation and deflation
- Radio opaque line and markings provided to facilitate identification of tube position
- Universal 15 mm connector at proximal end
- E.T.O Sterile and Pyrogen free. For single patient use

Packing

- Sterile, Individually packed in peelable pouch / blister packed.
- Box of 25; Master Carton of 300



Gross: 14.26kg (300 pcs)



660 x 380 x 500mm (5.0 to 10mm)

SIZES - CUFFED

5.0 mm, 5.5 mm, 6.0 mm, 6.5 mm, 7.0 mm, 7.5 mm, 8.0 mm, 8.5 mm, 9.0 mm, 9.5 mm, 10 mm